




MIXPAC™

T-Mixer™ Colibri™ Plus Breakable

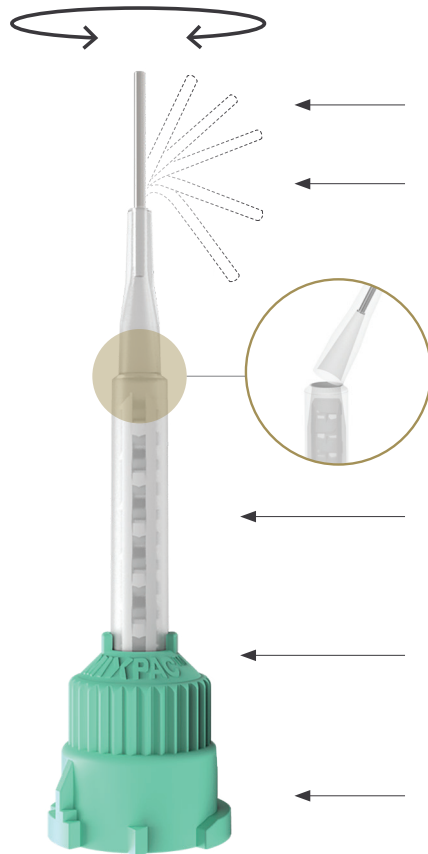

360° rotatable

Smooth rounded design

Bendable integrated cannula

Allows flexible and highly precise dispensing, even in difficult to access area.

Integrated breaking zone

For quick and easy application

T-Mixer technology

Up to 40% less material waste

MIXPAC Stamp

Genuine MIXPAC mixing tips have the MIXPAC brand stamped.

T-Mixer Colibri Plus Breakable


One color for all-in-one solution

Product Advantages

Integrated smoothly rounded stainless steel cannula	<ul style="list-style-type: none"> • Bendable and 360° rotatable - easy access to difficult areas. • Constant flow of material, even when the cannula is bent. • No additional intraoral tip is required. • Minimized risk of soft tissue trauma due to smooth, rounded cannula.
Predetermined breaking zone	<ul style="list-style-type: none"> • Front part of mixing tip is removable fast and easy. • Application of large amounts of material.
T-Mixer technology	<ul style="list-style-type: none"> • Up to 40% less material waste • Homogenous mixing results
All-in-one solution	<ul style="list-style-type: none"> • Impression of root canal • Impression of preparations and cavities • Cementation of root post • Core build-up



ALL-IN-ONE SOLUTION

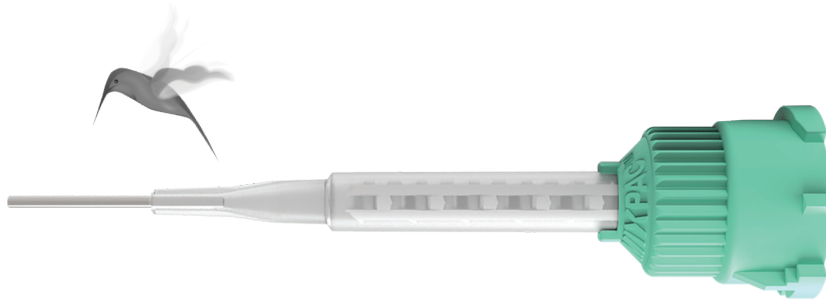

Ratio	Description	Cannula					Element ø (mm)	Length (mm)*	Waste Vol. (mL)	Order No.	Article
		MOQ Box	Outer ø (mm)	Inner ø (mm)	Gauge	Elements					
1:1, 2:1	Mint, natural elements	250	1.1	0.8	19	12	4.2	63.8+15	0.62	165788	MBT 4.2-12-N11B

* length mixing tip + length cannula

Your added value	
Offer your customers a complete system with the latest technology	<ul style="list-style-type: none"> • Fulfills user requirement toward simplification • All-in-one solution • Less inventory required • Be a pioneer of a new innovative mixing solution • Increased product brand awareness due to differentiation with a new mixing tip
A partner you can trust	<ul style="list-style-type: none"> • Customized bag packaging on request • Support in marketing and communication to end-users

What dentists say!



Dr. Eric Van Dooren
Antwerp, Belgium

"Being able to use almost the same mixing tip for the majority of clinical situations is an absolute must to minimize inventory and storage requirements in my dental practice. It's easier to use an all-in-one mixing system without the need to attach an intra oral tip in order to receive the key advantages of accurate application, while the anterior part of the mixing tip can be quickly and easily removed which allows for very fast dispensing/application."



Prof. Dr. Carlos Eduardo Francci
Sao Paulo, Brazil

"The T-Mixer Colibri Plus Breakable was designed to improve upon the existing innovative mixing tip technology for easier and faster delivery of the material with the integrated advantages of the T-Mixer technology."

medmix Switzerland AG | Ruetistrasse 7 | 9469 Haag, Switzerland
P +41 81 414 70 00 | www.medmix.swiss



Legal Notice: This brochure is a general presentation. The information contained in this publication is believed to be accurate and reliable but is not to be construed as implying any warranty or guarantee of performance. We waive any liability and indemnity for effects resulting from its application. Please, contact us for a description of the warranties and guarantees offered with our products. Directions for use and safety will be given separately. All information herein is subject to change without notice. MIXPAC, T-Mixer, and Colibri are trademarks of medmix Switzerland AG. medmix is a trademark of medmix AG. © medmix Switzerland AG 2023



DENTAL

DE-MX-DAT-OEM-Colibri Breakable 02.23 en