


# transcodent™ Dental Needles


German engineered products for comfort and safety

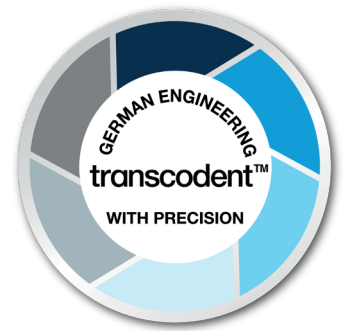

- Painless Steel™ Dental Injection Needles
- Luminject™ Dental Injection Needles
- Injection Needles
- Intraligamentalia™
- Suction Needles
- Endo Evacuation Needles
- Endo Irrigation Needles
- Irrigation Needles
- Paste Needles

## Cylinder ampoules

- Plastic hub, for use with cylinder ampoule syringes
- Siliconized cannula
- G 30 for the intraligamentary anaesthesia
- Precisely threaded hub for secure fit and easy handling
- EO sterilized

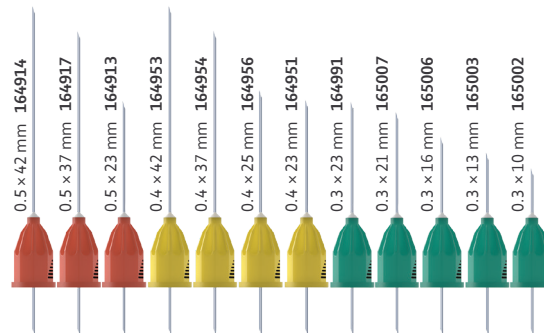

## Painless Steel™ Dental Injection Needles



**Cylinder ampoules**  
 Sizes: G 25 / 0.5 mm, G 27 / 0.4 mm, G 30 / 0.3 mm  
 Colors: red, yellow, green  
 Packaging: 100 pieces



Increased flow.  
Enlarged lumen

## Cylinder ampoules

- For intraligamentary and infiltration anaesthesia. Upper jaw-, lower jaw- and anterior teeth area
  - Increased flow due to enlarged lumen
  - Thin Wall
- Luer Lock**
- Luer Lock hub for single-use syringe

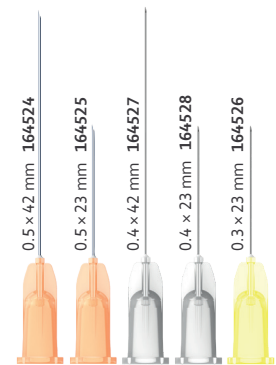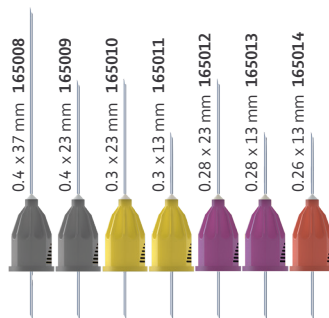

## Luminject™ Dental Injection Needles

**Cylinder ampoules Thin Wall**  
 Sizes: G 27 / 0.4 mm, G 30 / 0.3 mm, 0.26 mm, 0.28 mm  
 Colors: grey, yellow, violet, red  
 Packaging: 100 pieces

## Injection Needles

**Luer / Luer Lock**  
 Sizes: G 25 / 0.5 mm, G 27 / 0.4 mm, G 30 / 0.3 mm  
 Colors: red, grey, yellow  
 Packaging: 100 pieces


- Cylinder ampoule syringe for the intraligamentary and intraseptal anaesthesia
- High quality stainless steel autoclavable and sterilizable
- Multiple quadrants within one session possible
- Comfortable and precise dosage: ≈0.06 mL per click
- Pressure limitation



164604

## Intraligamentalia™

**Warranty:** 2 Years  
**Material:** Stainless Steel  
**Packaging:** 1 piece

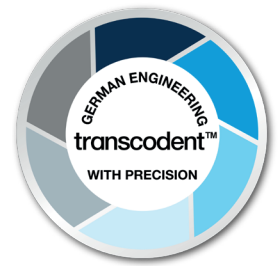

# transcodent™ Dental Needles



Annealed steel.  
Easy to bend

- Soft annealed cannula tube
- Easy to bend
- For suctions in places that are difficult to access, e.g. during endodontic or surgical procedures

\* TW (Thin Wall) cannula



Apical side vent



Laser marks

- Apical side vent
- For suctions near the apex
- 4 laser markings



## Suction Needles

### Luer Lock

Sizes: G 14 / 2.1mm, G 16 / 1.6mm, G 18 / 1.2 mm, G 20 / 0.9mm  
Colors: green, white, pink, yellow  
Packaging: 100 pieces

## Endo Evacuation Needles

### Luer Lock

Sizes: G 27 / 0.4mm, G 30 / 0.3mm  
Colors: yellow, grey  
Packaging: 25 pieces

- Irrigation of the full length of the root canal
- Luer Lock base, single packed
- Thin Wall – particularly thin-walled needle for more flexibility
- EO sterilized

## Double-sided vents

- The Original: Irrigation Needles with 2 lateral openings

- Penetration of the apex is reduced due to the closed, rounded tip

## Single-sided vent

- For gentle but effective irrigation of root canal. Penetration of the apex is avoided due to the lateral vent and rounded end



Double-sided vents



Single-sided vent



## Endo Irrigation Needles

### Luer Lock Thin Wall

Sizes: G 23 / 0.6mm, G 25 / 0.5mm, G 27 / 0.4mm, G 30 / 0.3mm, G 31 / 0.28mm  
Colors: blue, red, grey, yellow, pink  
Packaging: 25 pieces

### Luer Lock Thin Wall

Sizes: G 25 / 0.5mm / G 30 / 0.3mm  
Colors: red, yellow  
Packaging: 25 pieces

## Luer Lock

- Irrigation of periodontal pockets and root canal
- Suitable for Luer and Luer Lock syringes
- Suitable for single-use syringe
- EO sterilized



## Irrigation Needles

### Luer / Luer Lock

Sizes: G 25 / 0.5mm, G 27 / 0.4mm, G 30 / 0.3mm  
Colors: yellow, blue, green  
Packaging: 100 pieces

- Rounded tip: for the application in periodontal pockets

- EO sterilized

## Cylinder ampoules

- With plastic base
- Suitable for common cylinder ampoule syringes
- Easy to bend

## Luer Lock

- Suitable for Luer and Luer Lock syringes



### Cylinder ampoules

Sizes: G 20 / 0.9 mm, G 21 / 0.8 mm, G 23 / 0.6 mm  
Colors: red, green, blue  
Packaging: 100 pieces



### Luer / Luer Lock

Sizes: G 20 / 0.9 mm, G 23 / 0.6 mm  
Colors: red, yellow, blue  
Packaging: 100 pieces

## Paste Needles

Item designation: e.g. (0.5 x 42 mm 162071)  
(Diameter x length mm Article No.)

medmix Deutschland GmbH | Zum Kesselort 53 | 24149 Kiel | Germany www.medmix.swiss

Legal Notice: This brochure is a general presentation. The information contained in this publication is believed to be accurate and reliable but is not to be construed as implying any warranty or guarantee of performance. We waive any liability and indemnity for effects resulting from its application. Please, contact us for a description of the warranties and guarantees offered with our products. Directions for use and safety will be given separately. All information herein is subject to change without notice. medmix is a trademark of medmix AG, and transcodent, Painless Steel, Luminject and Intraligamentalia are trademarks of medmix Switzerland AG. Copyright © 2023 medmix Switzerland AG. All rights reserved.



DENTAL