

transodent™ Dental Needles

German engineered products for comfort and safety




- Painless Steel™ Dental Injection Needles
- Suction Needles
- Endo Evacuation Needles
- Flexible Steel™ Irrigation Needles



Cylinder ampoules

- Siliconized cannula with 3-edge lancet grinding provides smooth entry which minimizes pain and discomfort to the patient
- Hub bevel mark facilitates orientation of the bevel
- Threaded hub for easy and safe handling

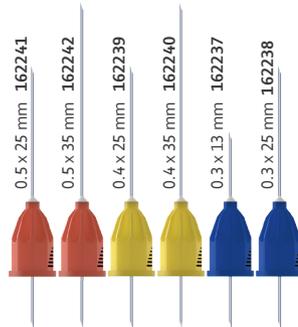

Painless Steel™ Dental Injection Needles	Cylinder ampoules Sizes: 25 G / 0.5 mm, 27 G / 0.4 mm, 30 G / 0.3 mm Colors: See picture Packaging: 100 Pieces
---	--



Annealed steel.
Easy to bend

◀39▶
mm

- Soft annealed cannula tube
- Easy to bend
- For suctions in places that are difficult to access, e.g. during endodontic or surgical procedures

* TW (Thin Wall) cannula



Apical side vent

◀32▶
mm



Laser Marks

- Apical side vent
- For suctions near the apex
- 4 laser markings

NEW



Suction Needles	Luer Lock Sizes: 14 G / 2.1mm, 16 G / 1.6mm, 18 G / 1.2mm, 20 G / 0.9mm Colors: See picture Packaging: 100 Pieces
------------------------	---

Endo Evacuation Needles	Luer Lock Sizes: 27 G / 0.4mm, 30 G / 0.3mm Colors: See picture Packaging: 25 Pieces
--------------------------------	--



Double-sided vents

◀25▶
mm

- Irrigation of the full length of the root canal
- Luer Lock base, single packed
- Thin Wall – particularly thin-walled needle for more flexibility
- Ethylene oxide sterilized
- The Original: Irrigation Needles with 2 lateral openings
- Penetration of the apex is reduced due to the closed, rounded tip



Flexible Steel™ Irrigation Needles	Endo Irrigation Luer Lock Thin Wall Sizes: 23 G / 0.6 mm, 25 G / 0.5 mm, 27 G / 0.4 mm, 30 G / 0.3 mm Colors: See picture Packaging: 40 Pieces
---	--



medmix US Inc., 8 Willow Street, Salem, New Hampshire 03079 | www.medmix.swiss/en-US

Legal Notice: This brochure is a general presentation. The information contained in this publication is believed to be accurate and reliable but is not to be construed as implying any warranty or guarantee of performance. We waive any liability and indemnity for effects resulting from its application. Please, contact us for a description of the warranties and guarantees offered with our products. Directions for use and safety will be given separately. All information herein is subject to change without notice. medmix is a trademark of medmix AG, and transodent, Painless Steel, Luminject and Intraligamentalia are trademarks of medmix Switzerland AG. Copyright © 2023 medmix Switzerland AG. All rights reserved.



DENTAL