


DENTAL

Product Overview

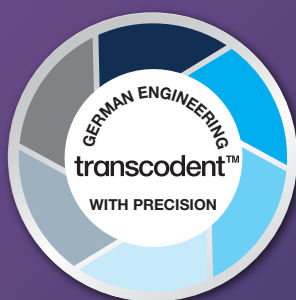

Imagine new and intelligent system solutions for mixing, dispensing and application of dental materials.

Welcome to medmix.




Swiss quality for innovative dentistry

Precision is crucial in dentistry, and dental professionals are constantly looking for new and intelligent solutions to reliably achieve the most accurate results possible. As a global leader in high-precision application systems, medmix provides its customers in the dental industry with uncompromising quality and excellent service in the form of the MIXPAC™ and transcotent™ brands.



MIXPAC is the original for innovative one- and two-component system solutions used to accurately mix, dispense and apply different types of dental materials for a variety of applications.

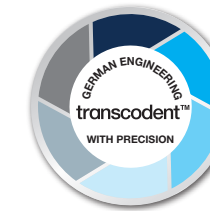

transcotent offers a high-quality range of sterile injection needles for anesthesia as well as needles for endo irrigation and evacuation applications.

Facts + figures


Headquarters:	Baar, CH
Production sites:	13 worldwide
Employees:	Over 2,000
Patents:	Over 900 active

Painless Injection
Effective Irrigation
Easy Evacuation
the
Original

The transcodent™ brand

Our German quality products provides a wide range of sterile, state-of-the-art needles for dental injection, endo irrigation and suction applications.

We offer highly specialized needles designed for patient comfort, better and more painless dental procedures and safety.

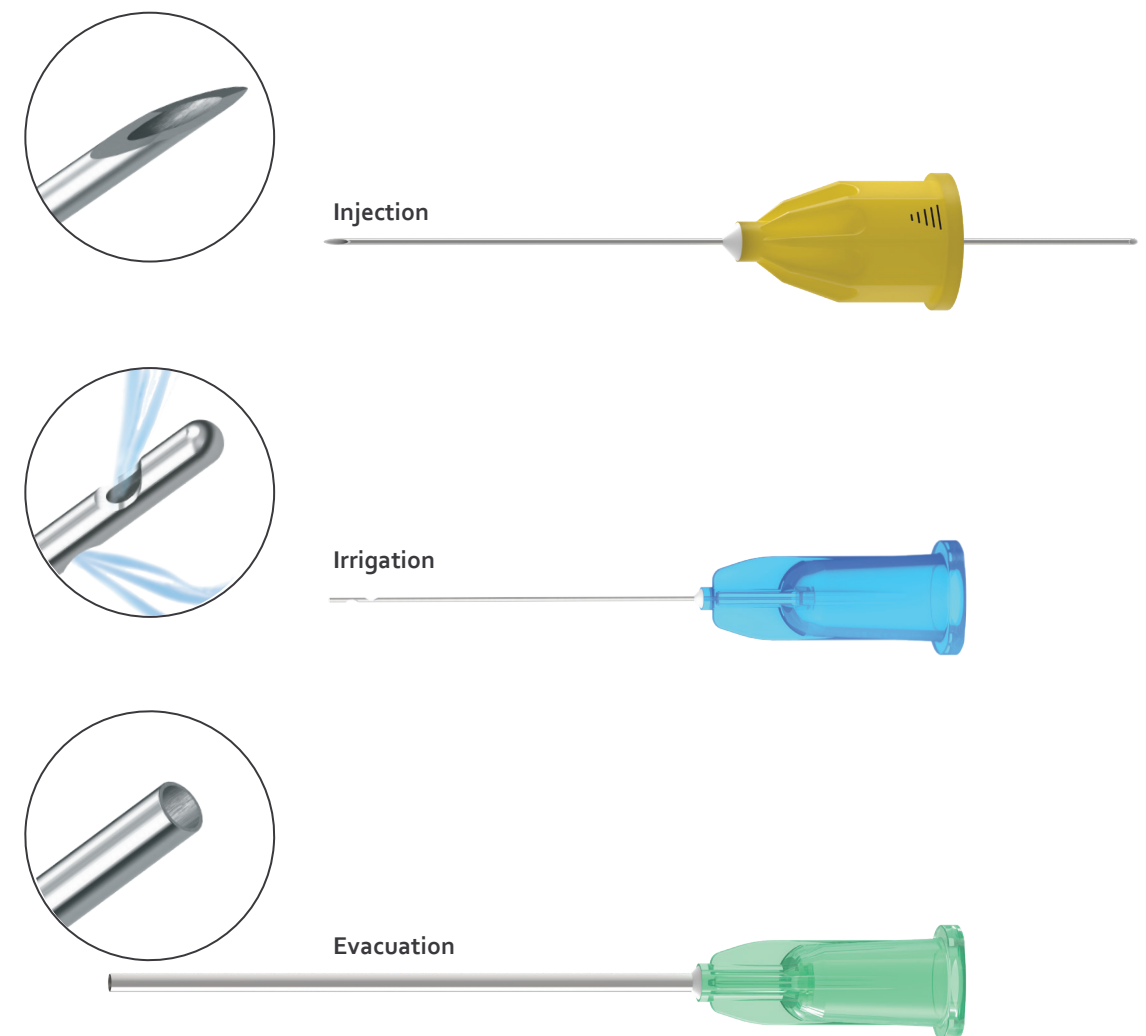


Dental Products

Dental Injection Needles	9
Painless Steel™	10
Irrigation Needles	12
Flexible Steel™	14
Suction Needles	16
Suction Needles	18
Endo Evacuation Needles	20
Barriers	22
Protective Sleeves 50 mL	23
Systems 5 mL	
Systems 1 mL	

Highly specialized dental products for patient comfort and more painless dental procedures and safety.

transcodent™

Injection Needles

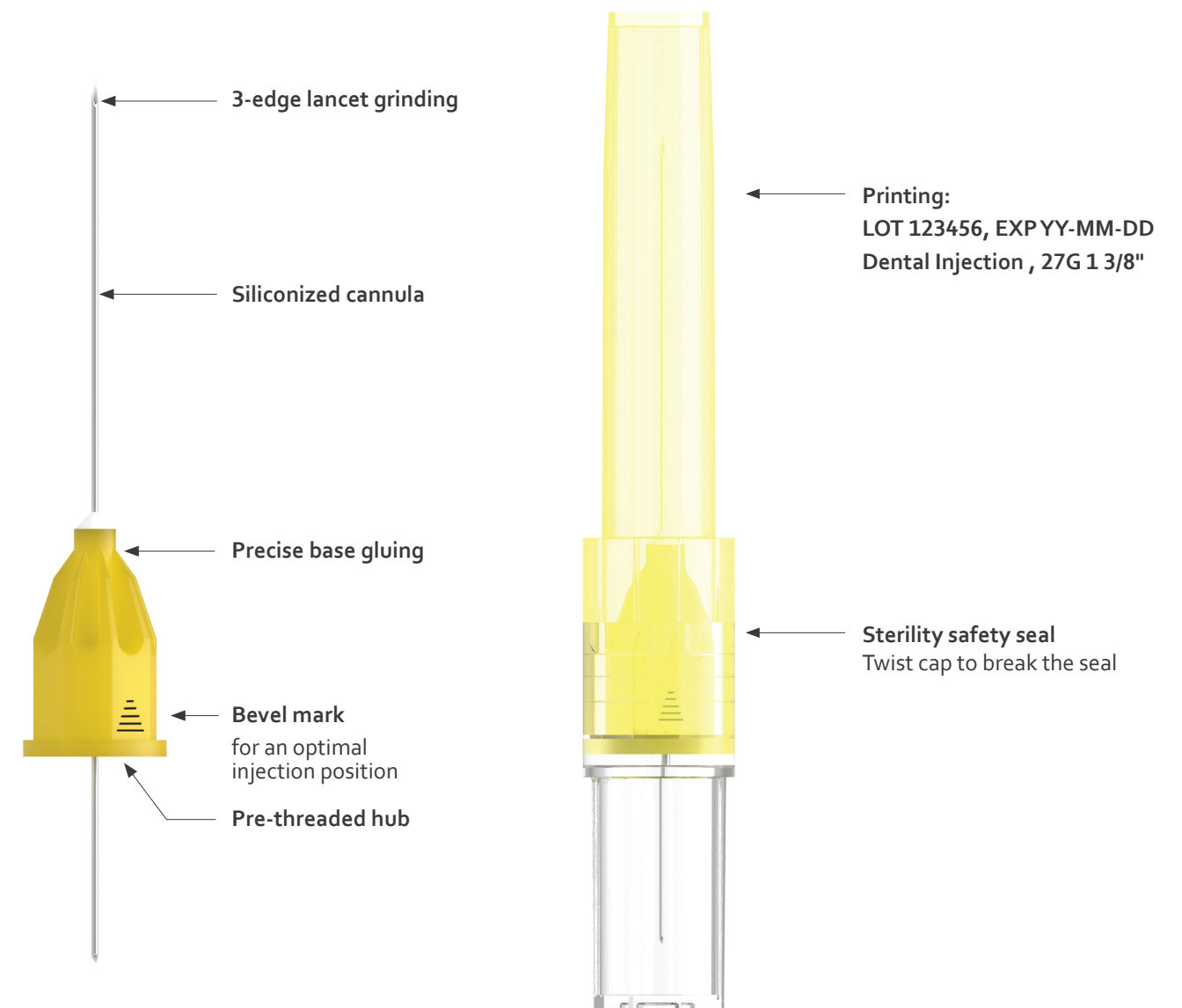


FEATURES

- 3-edge-lancet-grinding
- Siliconized cannula
- Bevel mark
- Threaded hub
- Precise base gluing
- Sterilized with ethylene-oxide
- Heat-sealed needle-container
- Suitable for cylinder ampoule syringes
- Precisely grinded cannula-back-end

Painless Steel™

DENTAL INJECTION NEEDLES

transcodent™ Painless Steel Dental Injection Needles are designed for effective and safe delivery of anesthetic solution to ensure optimal patient comfort.

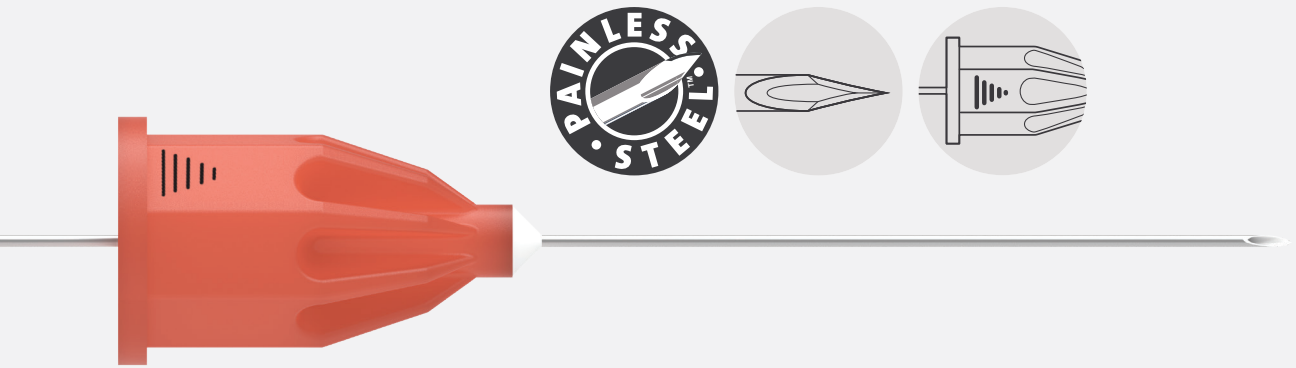

PRODUCT FEATURES AND BENEFITS

- **Siliconized cannula with 3-edge lancet grinding** provides smooth entry which minimizes pain and discomfort to the patient
- **Hub bevel mark** facilitates orientation of the bevel
- **Threaded hub** for easy and safe handling

Painless Steel™


Imperial



Order No.	Description	Size	Color	Qty/Box
162238	30 G Short	0.3 x 25 mm	blue	100
162237	30 G Extra Short	0.3 x 13 mm	blue	100
162240	27 G Long	0.4 x 35 mm	yellow	100
162239	27 G Short	0.4 x 25 mm	yellow	100
162242	25 G Long	0.5 x 35 mm	red	100
162241	25 G Short	0.5 x 25 mm	red	100

Endo Irrigation Needles are irrigation needles with double-sided vents for gentle but effective irrigation of root canals.

Irrigation Needles

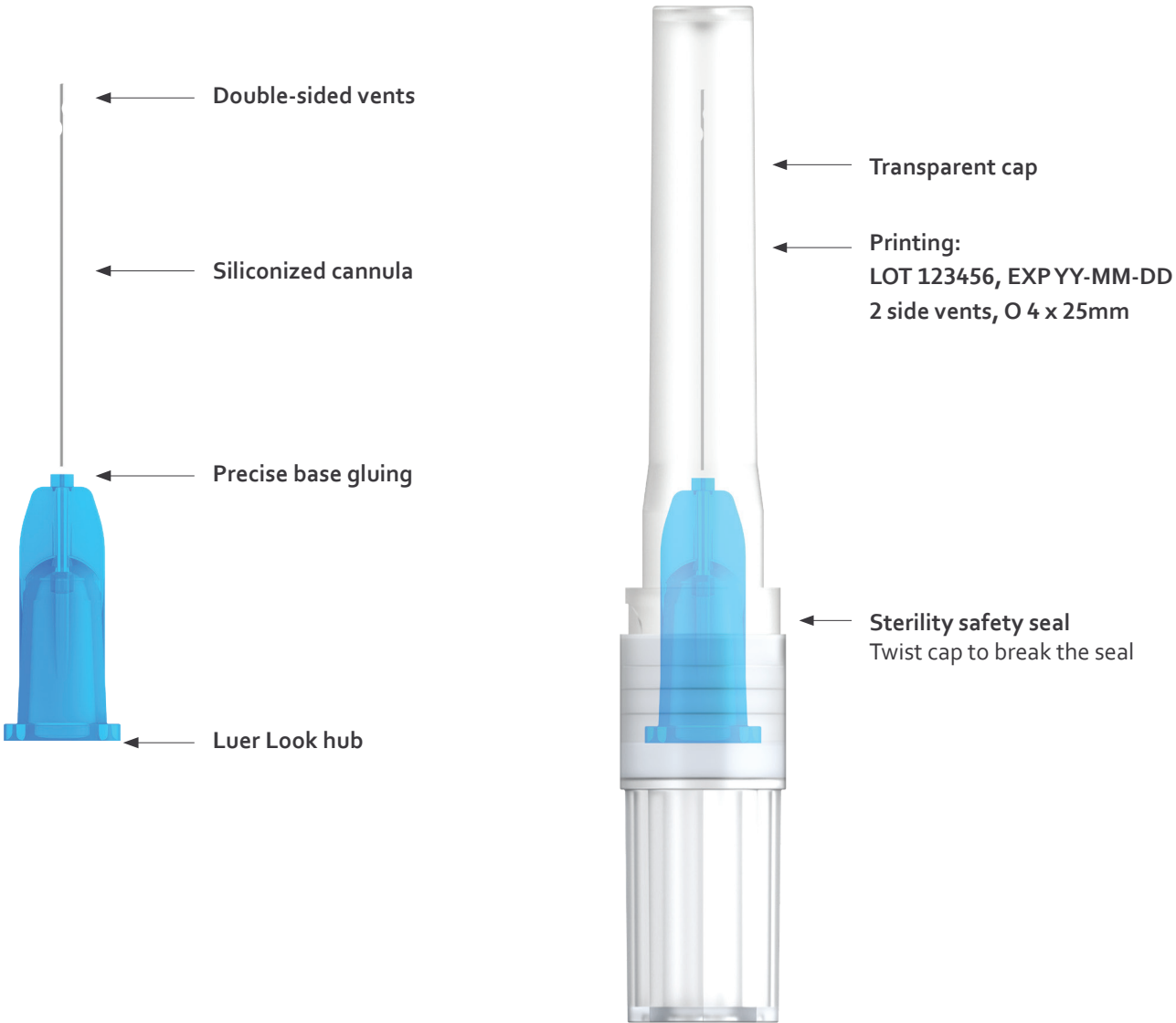


FEATURES

- Double-sided vents
- Sterilized with ethylene-oxide
- Easy to bend for perfect access to the cavity
- Luer Lock hub

Flexible Steel™


IRRIGATION NEEDLES

transcodent™ Flexible Steel Irrigation Needles are double-sided vent irrigation needles with Luer Lock that increase safety and efficacy of endodontic irrigation. The full working length can be easier achieved and effectively irrigated even in curved root canals.

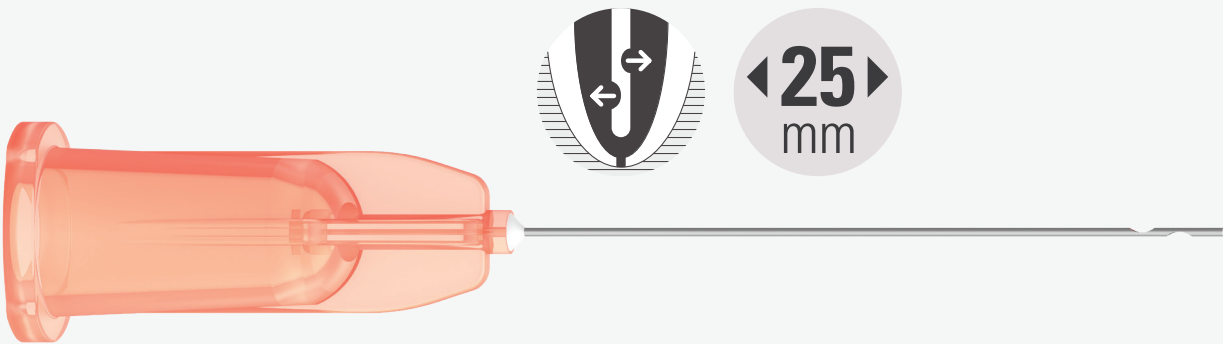

PRODUCT FEATURES AND BENEFITS

- **Closed-rounded tip with double-sided vents** minimizes the risk of needle binding and wedging in the root canals, trauma to soft/hard tissue and extrusion of the irrigation solution beyond apical foramen
- **Thin-walled cannula** provides improved flow rate for optimal flushing effect
- **Bendable cannula** for increased flexibility and improved access

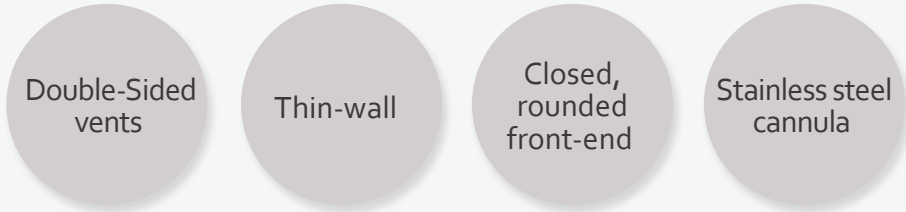

Flexible Steel™


Double-Sided Vents / Luer Lock



Order No.	Description	Size	Color	Qty/Box
162066	23G	0.6 x 25 mm	blue	40
162067	25G	0.5 x 25 mm	orange	40
162068	27G	0.4 x 25 mm	grey	40
162069	30G	0.3 x 25 mm	yellow	40

Designed for effective suction in areas with difficult access such as root canal or in a surgical site during oral surgery.

Suction Needles

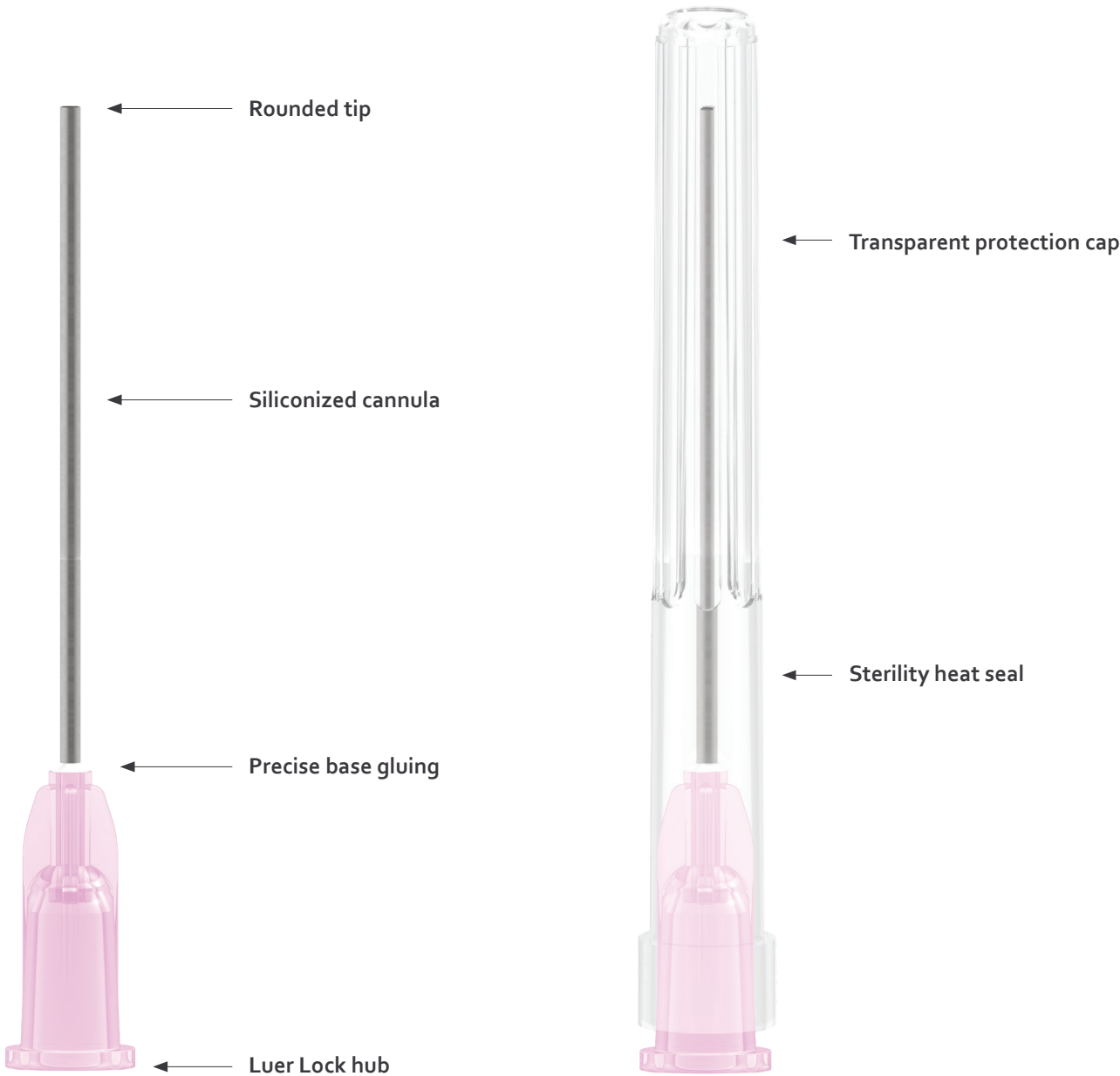



FEATURES

- Various gauge options for effective suctioning
- Laser markings to minimize apex trauma
- Sterilized with ethylene-oxide
- Easy to bend for perfect access to the cavity
- Suitable for cylinder ampoule syringes
- Luer Lock hub

Suction Needles

WITH LUER LOCK

transcodent™ Suction Needles thread are designed for effective suction in areas with difficult access such as root canal or in a surgical site during oral surgery. The annealed cannula is easy to bend.

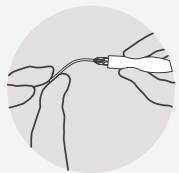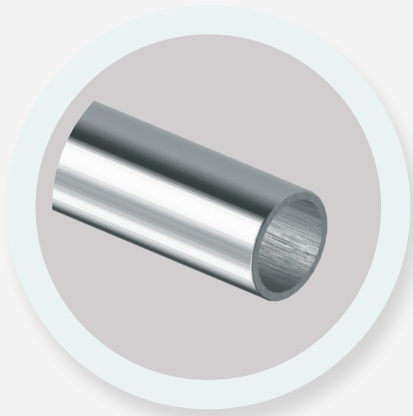


◀ 39 ▶
mm



PRODUCT FEATURES AND BENEFITS

- **Bendable, Thin-wall, Thin annealed steel cannula** enhances field of view and accessibility
- **Rounded tip** minimizes trauma to soft/hard tissue
- For use with the **Luer Lock Handpieces**.

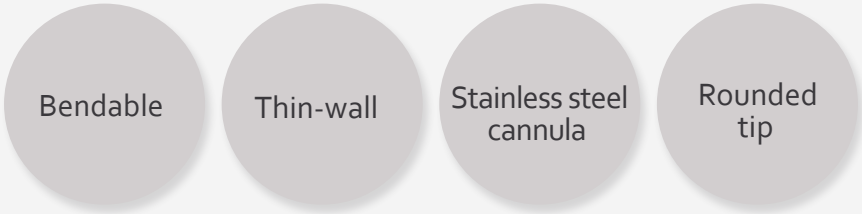

Suction Needles

Luer Lock




Order No.	Description	Size	Color	Qty/Box
162062	14G	2.1 x 39 mm	green	100
162063	16G	1.6 x 39 mm	white	100
162064	18G	1.2 x 39 mm	pink	100
162065	20G	0.9 x 39 mm	yellow	100

Endo Evacuation Needles


WITH LUER LOCK

transcodent™ Endo Evacuation Needles are designed to reduce trauma to the apex while providing effective suctioning in difficult access areas. The thin wall needle is remarkably flexible, ensuring efficient irrigation even in curved root canals.

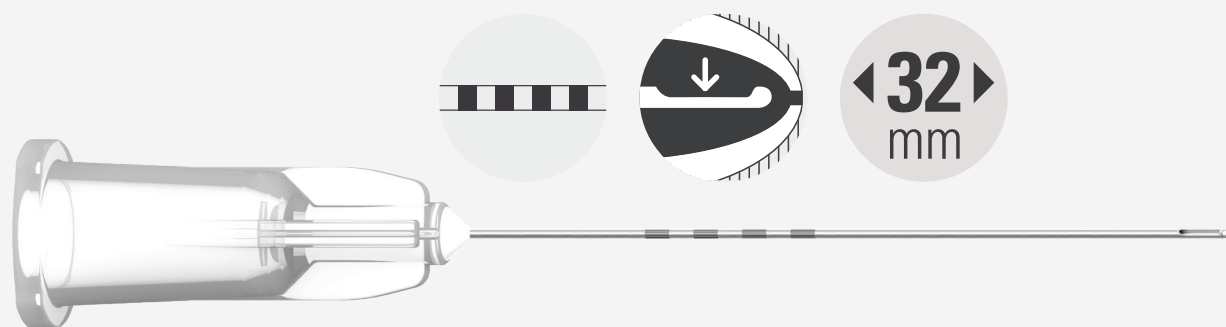

PRODUCT FEATURES AND BENEFITS

- **Apical side vent** for suctioning near the apex
- **4 laser markings** to reduce trauma to the apex
- **Luer lock hub** for use with the Luer Lock Handpieces

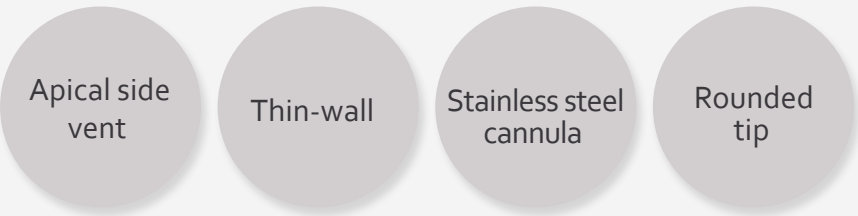


Endo Evacuation Needles

Luer Lock



Order No.	Description	Size	Color	Qty/Box
162054	27 G, sterile	0.4 x 32 mm	grey	25
162055	30 G, sterile	0.3 x 32 mm	yellow	25

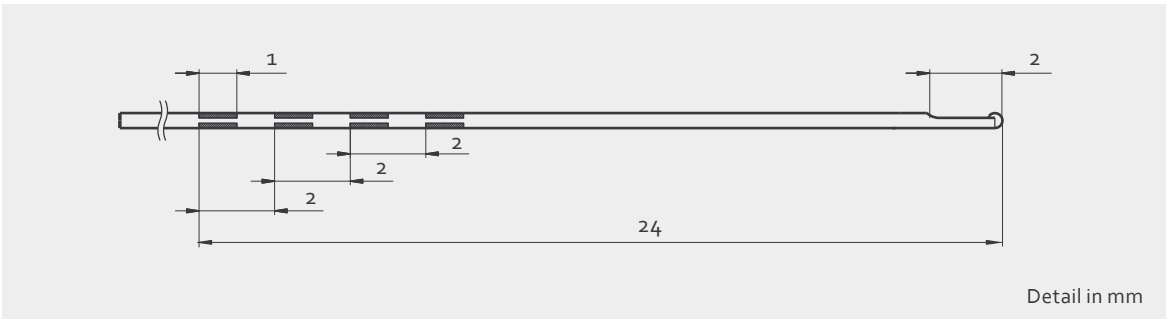

Barriers

FOR MIXPAC™ APPLICATION SYSTEMS



- Disposable single-use barrier
- Customized ergonomic design
- Tried and tested material
- Intuitive and user-friendly

MIXPAC Protective Sleeves provide convenient and effective barrier protection of multi-use application systems to protect the dentist, staff and patient against infectious disease during clinical use.

MIXPAC™

Protective Sleeves

Available for 50 mL, 5 mL, 1 mL MIXPAC cartridge systems



- 50 mL Systems
- 5 mL Systems
- 1 mL Systems

Order No.	Description	Material	Color	L (mm)	W (mm)	Qty/Box
168058	50 mL	PE-LD	Transparent	150	68	300
168057	5 mL	PE-LD	Transparent	220	51	300
168056	1 mL	PE-LD	Transparent	152	34	300

PE-LD = Low Density Polyethylene

Manufacturer: **Minitube AB** / Byvägen 44
835 96 Trångsviken / Sweden

FOR MORE INFORMATION, PLEASE CONTACT YOUR DEALER.

USA

medmix US Inc.
8 Willow Street
Salem, NH 03079, USA
P +1 603 893 2727
AME_Dental_Sales@medmix.com

www.medmix.swiss
www.medmix.swiss/en-US

EUROPE

medmix Switzerland AG
Ruetistrasse 7
9469 Haag, Switzerland
P +41 81 414 70 00

medmix Deutschland GmbH
Zum Kesselort 53
24149 Kiel, Deutschland
P +431 979 9790

ASIA

medmix China Ltd.
68 Bei Dou Lu, Minhang Shanghai
200245, China
P +86 21 6430 68 68



Legal Notice: This brochure is a general presentation. The information contained in this publication is believed to be accurate and reliable but is not to be construed as implying any warranty or guarantee of performance. We waive any liability and indemnity for effects resulting from its application. Please, contact us for a description of the warranties and guarantees offered with our products. Directions for use and safety will be given separately. All information herein is subject to change without notice. medmix is a trademark of medmix AG, and transcendent, Flexible Steel and Painless Steel, are trademarks of medmix Switzerland AG. Copyright © 2023 medmix Switzerland AG. All rights reserved.