


transcodent™ Dental Needles

German engineered products for comfort and safety




- Painless Steel™ Dental Injection Needles
- Suction Needles
- Endo Evacuation Needles
- Flexible Steel™ Irrigation Needles

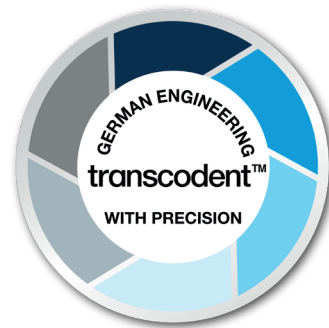

Cylinder ampoules

- Siliconized cannula with 3-edge lancet grinding provides smooth entry which minimizes pain and discomfort to the patient
- Hub bevel mark facilitates orientation of the bevel
- Threaded hub for easy and safe handling

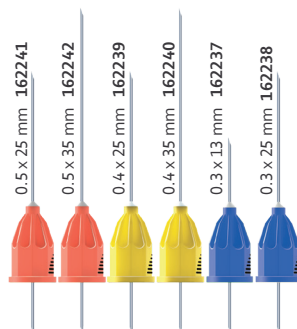

Painless Steel™ Dental Injection Needles

Cylinder ampoules

Sizes: 25 G / 0.5 mm, 27 G / 0.4 mm, 30 G / 0.3 mm
Colors: See picture
Packaging: 100 Pieces




Annealed steel.
Easy to bend

◀39▶
mm

- Soft annealed cannula tube
- Easy to bend
- For suctions in places that are difficult to access, e.g. during endodontic or surgical procedures

* TW (Thin Wall) cannula

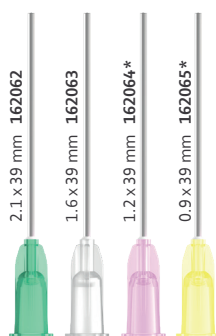

Apical side vent



Laser Marks

- Apical side vent for suctioning near the apex
- 4 laser markings to reduce trauma to the apex
- Luer lock hub

◀32▶
mm



NEW

Suction Needles

Luer Lock

Sizes: 14 G / 2.1mm, 16 G / 1.6mm, 18 G / 1.2mm, 20 G / 0.9mm
Colors: See picture
Packaging: 100 Pieces

Endo Evacuation Needles

Luer Lock


Sizes: 27 G / 0.4mm, 30 G / 0.3mm
Colors: See picture
Packaging: 25 Pieces



Double-sided vents

◀25▶
mm

- Irrigation of the full length of the root canal
- Luer Lock base, single packed
- Thin Wall – particularly thin-walled needle for more flexibility
- Ethylene oxide sterilized
- The Original: Irrigation Needles with 2 lateral openings
- Penetration of the apex is reduced due to the closed, rounded tip

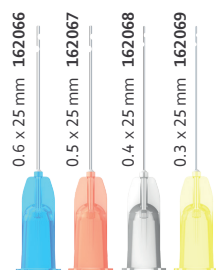

Flexible Steel™ Irrigation Needles

Endo Irrigation Luer Lock Thin Wall

Sizes: 23 G / 0.6 mm, 25 G / 0.5 mm, 27 G / 0.4 mm, 30 G / 0.3 mm
Colors: See picture
Packaging: 40 Pieces

medmix US Inc. | 8 Willow Street | Salem, NH 03079, USA | www.medmix.swiss/en-US

medmix Deutschland GmbH | Zum Kesselort 53 | 24149 Kiel, Deutschland

Legal Notice: This brochure is a general presentation. The information contained in this publication is believed to be accurate and reliable but is not to be construed as implying any warranty or guarantee of performance. We waive any liability and indemnity for effects resulting from its application. Please, contact us for a description of the warranties and guarantees offered with our products. Directions for use and safety will be given separately. All information herein is subject to change without notice. medmix is a trademark of medmix AG, and transcodent, Painless Steel, Luminject and Intraligamentaria are trademarks of medmix Switzerland AG. Copyright © 2024 medmix Switzerland AG. All rights reserved.



DENTAL