




**MIXPAC™**

# T-Mixer™ Colibri™ Plus Breakable

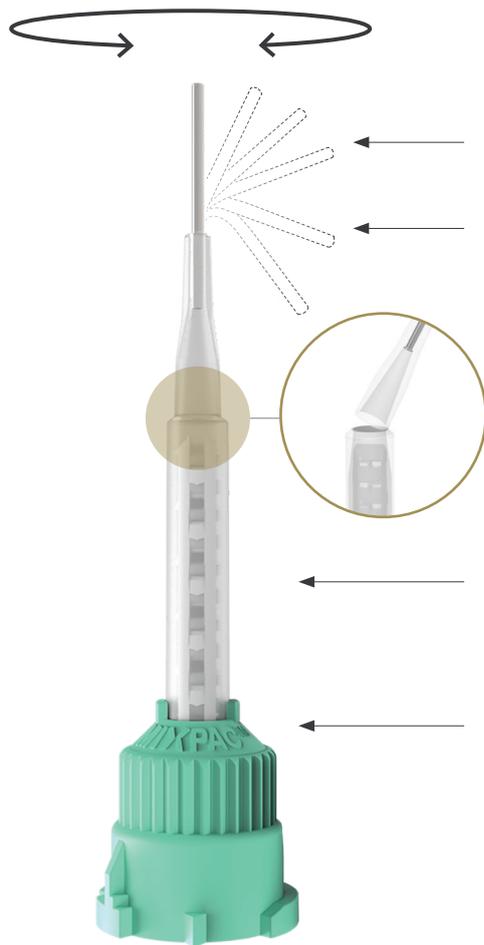

**360° rotatable cannula**  
for precise and customized application

**Smooth rounded design**  
minimizes trauma to soft tissue

**Bendable integrated cannula**  
allows flexible and highly precise dispensing,  
even in difficult-to-reach areas

**Integrated breaking zone**  
facilitates quick and easy application of large  
amounts of material

**T-Mixer technology**  
Up to 40% less material waste


**MIXPAC Stamp**  
genuine MIXPAC mixing tips have the MIXPAC  
brand stamped on the top edge

## Product Advantages

Integrated smoothly rounded stainless steel cannula	<ul style="list-style-type: none"> <li>• Bendable and 360° rotatable - easy access to difficult-to-reach areas.</li> <li>• Constant flow of material, even when the cannula is bent.</li> <li>• No additional intraoral tip is required.</li> <li>• Minimized risk of soft tissue trauma due to smooth, rounded cannula.</li> </ul>
Predetermined breaking zone	<ul style="list-style-type: none"> <li>• Front part of mixing tip is removable fast and easy.</li> <li>• Application of large amounts of material.</li> </ul>
T-Mixer technology	<ul style="list-style-type: none"> <li>• Up to 40% less material waste.</li> <li>• Homogenous mixing results.</li> </ul>
All-in-one solution	<ul style="list-style-type: none"> <li>• Impression of root canal.</li> <li>• Impression of preparations.</li> <li>• Cementation of root post.</li> <li>• Core build-ups.</li> </ul>



# ALL-IN-ONE SOLUTION

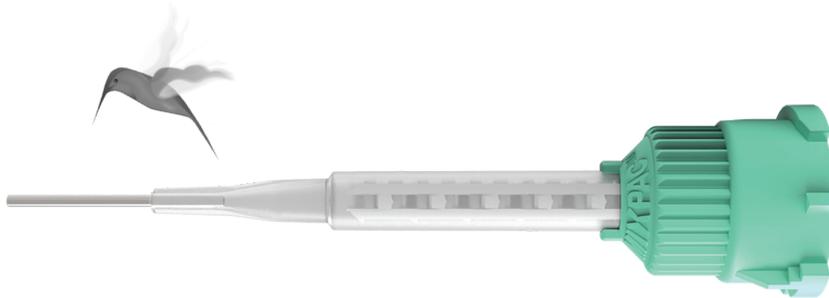

Ratio	Description	Cannula				Element ø (mm)	Length (mm)*	Waste Vol. (mL)	Order No.	Article
		MOQ Box	Outer ø (mm)	Inner ø (mm)	Gauge Elements					
1:1, 2:1	Mint, natural elements	250	1.1	0.8	19 12	4.2	63.8+15	0.62	165788	MBT 4.2-12-N11B

\* length mixing tip + length cannula

## Your added value

Offer your customers a complete system with the latest technology

- Fulfills user requirement toward simplification.
- All-in-one solution.
- Less inventory required.
- Be a pioneer of a new innovative mixing solution.
- Increased product brand awareness due to differentiation with a new mixing tip.

A partner you can trust

- Customized bag packaging on request.
- Support in marketing and communication to end-users.

## What dentists say!



**Dr. Eric Van Dooren**  
Antwerp, Belgium

*"Being able to use almost the same mixing tip for the majority of clinical situations is an absolute must to minimize inventory and storage requirements in my dental practice. It's easier to use an all-in-one mixing system without the need to attach an intra oral tip in order to receive the key advantages of accurate application, while the anterior part of the mixing tip can be quickly and easily removed which allows for very fast dispensing/application."*



**Prof. Dr. Carlos Eduardo Francci**  
Sao Paulo, Brazil

*"The T-Mixer Colibri Plus Breakable was designed to improve upon the existing innovative mixing tip technology for easier and faster delivery of the material with the integrated advantages of the T-Mixer technology."*

medmix Switzerland AG | Ruetistrasse 7 | 9469 Haag, Switzerland  
P +41 81 414 70 00 | [www.medmix.swiss](http://www.medmix.swiss)



Legal Notice: This brochure is a general presentation. The information contained in this publication is believed to be accurate and reliable but is not to be construed as implying any warranty or guarantee of performance. We waive any liability and indemnity for effects resulting from its application. Please, contact us for a description of the warranties and guarantees offered with our products. Directions for use and safety will be given separately. All information herein is subject to change without notice. MIXPAC, T-Mixer, and Colibri are trademarks of medmix Switzerland AG. medmix is a trademark of medmix AG. © medmix Switzerland AG 2023



DENTAL